**VIROMED**

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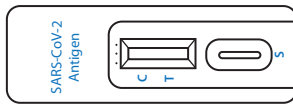
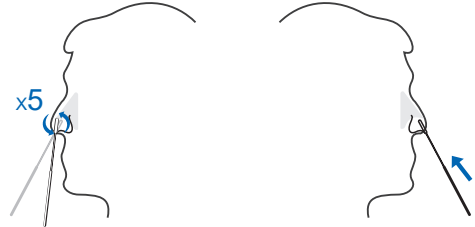
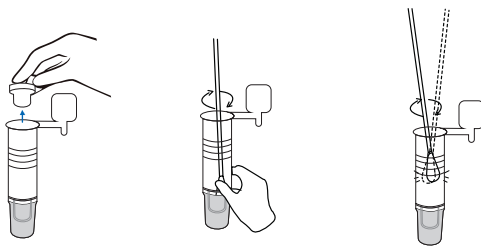
STEP-BY-STEP-INSTRUCTION

1

Open the sealed pouch and remove the test cassette. Lay it face up on a clean, dry and flat surface.

2

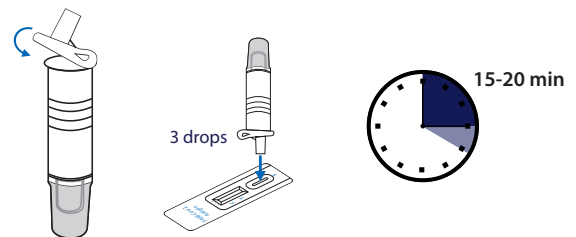
Gently, insert the entire absorbent tip of the swab (around 1.5 cm) into your nostril. Slowly, rotate the swab in a circular against the inside walls of your nostril 5 times or more. Be sure to collect any nasal drainage that maybe present on the swab. Gently remove the swab. Use the same swab to repeat steps in the other nostril and slowly, take out the swab.

3

Take sample tube with prefilled sample extraction buffer and remove the white cap of the sample tube.

Insert the swab into the sample tube prefilled with extraction buffer. Mix well and squeeze the swab 10-15 times by compressing the walls of the tube against the swab.

Roll the swab head against the inner wall of the tubes as you remove it. Try to release as much liquid as possible. Dispose of the used swab in accordance with your local regulations.

4

Close the cap of the sample tube.

Add 3 full drops of the mixed solution vertically into the sample well (S) of the test cassette.

Read the result 15-20 minutes after adding the sample. Result got after 20 minutes is invalid.

5

Positive

Two colored bands appear on the membrane. One band appears in the control region (C) and another band appears in the test region (T).

The test result means that SARS-CoV-2 antigen are detectable in your sample. The detection of these antigens indicates with a high probability of infection with the novel coronavirus.

Please stay at home and contact your physician or the responsible health authority immediately to get information on how to proceed further.

***Note:** The thickness of the line is insignificant; any reddish color in the Test line (T) should be considered a positive result. The positive test result must be confirmed by PCR.

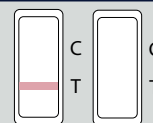
Negative

Only one colored band appears, in the control region (C). No apparent colored band appears in the test region (T).

The test result indicates that there is no or too little SARS-CoV-2 Antigen in the sample and at the current time there is probably no infection with the novel coronavirus.

A negative result does not preclude SARS-CoV-2 infection, so please stay at home if you have clinical symptoms or if you have a well-founded suspicion and contact physician or responsible health authority to get information on how to proceed further.

False negative results can be from incorrect sampling, incorrect execution of the test, or insufficient virus in the sample.

Invalid

If there is no Control line (C) or only a Test line (T) in the result window, the test did not run correctly and the results are not valid.

It is important that you carefully follow the instructions for the test. You should test again with a new sample and a new test.